

G340

Aligner box with a lid that snaps shut. Available in various 2 colour combinations. 2 aligners can fit into it easily. Measurements 10x7,6 h 2 cm.




G342

Aligner box with a lid that snaps shut. Available in various 2 colour combinations, with a place for the chewie and extraction hook. Measurements 10x7,6 h 2 cm.




G345

Aligner box with a lid that snaps shut. Available in various 2 colour combinations, with a mirror and a place for the chewie and extraction hook. Measurements 10x7,6 h 2 cm.

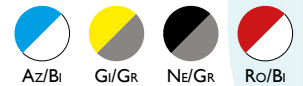

G43

Small retainer box (for bite) with ventilation holes various colours, with hole for cord. Available in 10 different colours. 8,4x8,3 cm, h 2,5 cm.

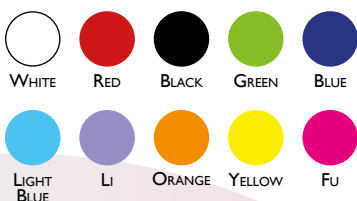

MEL-TRAY

Tray made of melamine, 28x18 cm, with molded spaces for instruments, trays and mixing wells sterilizable at 121°.



The MELAMINE Melochem we use complies with the compositional requirements needed by the Regulation 10/2011 / EC concerning materials intended to come into contact with food. The material is durable, robust, resistant to solvents and acids and sterilizable, therefore reusable.

Melamine resembles plastic but it is not plastic: it can be defined as a sustainable material, as it is a BPA-free thermosetting resin, i.e. free of bisphenol, and does not derive from hydrocarbons.

It is not porous, it inhibits the growth and presence of bacteria, viruses and fungal substances and is non deformable. Usually used in the food sector with products such as plates, cups and glasses. It is now also suitable for sanitary use with instrument trays and other devices in the dental sector. It can be sterilized at 121°C or with cycle B at 135°C. It does not break, does not deteriorate and does not deform.

MEL-52

Medium melamine tray container cm 20.7x14.4, perfect support for disposable tray code C52.



MEL-15

Small melamine tray cm 24x16.2 (internal base cm 19.5x12.5). Ideal container for dental instruments or support for disposable trays.



MEL-16

Tray made of melamine cm 28,4x18,3x1,7.



MEL-16.1

Melamine cover for tray MEL-16.

MEL-14

Small melamine tray cm 23x12 (internal base cm 19.5x9.5). Ideal container for dental instruments or support for disposable trays.



CANNULA "ZERO"

The "Zero" saliva ejector, which is called this because when used it **eliminates over 95% of the generated aerosol**, thanks to the flow rate of over **300 litres/minute**, it is made up of a **double mirror, double suction and selector**.

The selector allows you to select the side from which you prefer to work (without subjecting your wrist and elbow to inadequate movements and positions) and position it correctly in the mouth so as to divert the suction flow from one of the two sides of the saliva ejector and to close the opposite one to avoid that, once positioned, it could create a suction effect with many inconveniences for the patient. Turning the selector to "O" (Open) opens the suction side from which you choose to work. The suction side "C" (Closed) which

comes into contact with the mucous membrane or tongue is automatically excluded.


The material used is biocompatible, autoclavable, hypoallergenic, antitoxic and, above all, certified. As for the accessories (flexible tube and connectors) they can be useful to improve movements during the use of the saliva ejector itself, for this reason we have created light and space-saving connectors to help easy movement.

The tube is easy to get, there are different types, and just select a light one (such as those used for oxygen therapy).

It can be a useful support for those that useGBT therapy, scalign, dentists, preparations or surgery, etc. Available in 2 colours, black or white.

C582 Mirror diameter 22 mm, package with 1 Cannula "Zero" complete + 2 spare mirrors

C584 Mirror diameter 24 mm, package with 1 Cannula "Zero" complete + 2 spare mirrors

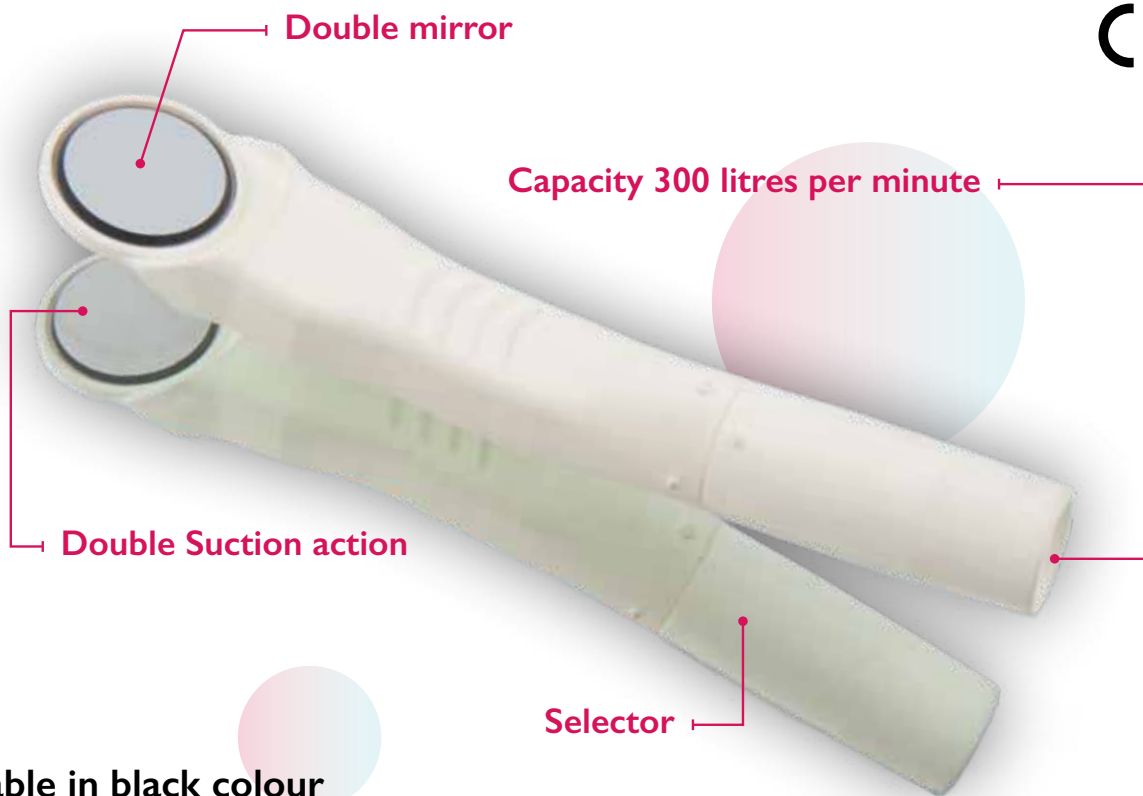

Available in black colour

WATCH
DEMO HERE



IMPORTANT

When changing the viewing side of the mirror it is necessary to select the suction side by acting on the selector. Then you have to turn the selector in the direction where you want to open the cannula suction

L42

Silicone autoclavable bur holder with 35 universal holes suitable for every kind of bur. Supplied with transparent, plastic cover. Available in 3 colours.



NANOSINT BRUSHES RANGE

Nanosint Microlayer brushes are specially made for the new techniques of micro layering and colouring on zirconium teeth to be finished. They are produced in Italy with new generation very thin synthetic bristles, which make them very similar to the kolinsky marten.

The main advantage is the resistance of the tip, which does not curl due to wear and tear and keeps its quality for a long time. Regarding the 3 larger models, PL-4 PL-9 PL-12, the reversible aluminum handle lets the soaked tip be put inside in order to maintain the correct humidity and be used as a cover



PL-00

Ultra thin for interdental and proximal nano staining and the handle in red wood - bristles length is 6 mm. Base thickness is 0.3 mm.



PL-0

Extra thin for micro colouring with a red wooden handle - bristle length is 8 mm. Base thickness is 0.6 mm.



PL-1

Thin to colour with a red wooden handle - bristle length is 10 mm. Base thickness is 1 mm.



PL-KIT

KIT of 6 Nano sint Microlayer brushes.

PL-4

Medium for micro layering and colouring, reversible aluminum handle - bristle length is 12 mm. Thickness at the base is 2.5 mm.



PL-9

Medium-large compact form for dense masses or traditional ceramic with reversible aluminum handle - bristle length is 13 mm. Base thickness is 4 mm.



PL-12

Big compact form for dense masses and traditional ceramic with a reversible aluminum handle - bristle length is 16 mm. Base thickness is 5 mm.



Larident Srl

Via Lamaneigra 12R (Ponte Rimassa) - I6030 Tribogna (Genova)
Tel. 0185 938068 - www.larident.it

